

JOHNLEN

MULTIWRAP MOIST & DRY HEATING PAD

INSTRUCTIONS



遠紅外線

多用途乾濕兩用電子熱敷墊

居家保健用 說明書

各部名稱

PARTS NAMES

軟性外套，可分離清洗
SOFT COVER, WITH SNAPS, IS
REMOVABLE & WASHABLE.

獨特拉鍊設計，可依外套不同的變化，
使用在人體各個部位
UNIQUE ZIPPER DESIGN FORMS ENVELOPE, SHAPING
PAD TO FIT VARIOUS PARTS OF YOUR BODY.

鬆緊帶定位適當位置
ELASTIC BANDS HOLD SOOTHING HEAT IN PLACE.

兩段控溫附指示燈
LEVER CONTROL WITH 2 HEAT
SETTINGS & ON INDICATOR
LIGHT.

黏扣帶
SELF-ADHERING HOOK & LOOP TAPE.

吸水海棉片可作濕式熱敷
ABSORBENT SPONGE PROVIDES
PENETRATING MOIST HEAT.

特點：

- 遠紅外線（波長 8 - 16 MICRON）滲透原理設計，具有促進血液循環，有益人體皮膚細胞之新陳代謝。
- 裝置恆溫器作溫度自動控制。
HI: 45°C ~ 50°C LO: 35°C ~ 40°C

效能：

- 酸痛之熱敷治療，如關節炎、風濕症、麻痺復健、神經痛、運動傷害等。
- 冬天時，手脚冰冷之保溫，幫助安眠。

CHARACTERISTIC：

- THE PENETRATING HEAT IS DESIGNED TO PROMOTE BLOOD CIRCULATION, HELPFUL METABOLISM OF HUMAN CELLS.
- INSTALL THERMOSTAT FOR HEAT SETTING.
HI: 45°C ~ 50°C LO: 35°C ~ 40°C

FUNCTION：

- WRAPS SOOTHING HEAT AROUND SORE OR STIFF AREA AFTER ACTIVE SPORTS & OTHER STRENUOUS ACTIVITIES, ARTHRITIS, RHEUMATISM, PARALYSIS HEALTH CARE NEURALGIA, ETC.
- KEEP HAND & FOOT WARM AT COLD WEATHER, HELPFUL SLEEP.

使用時注意事項

DANGER

- 不當使用會造成皮膚灼傷，敬請務必詳閱說明書。
- 睡覺時，切勿使用。
- 不能使用在小孩身上。
- 血液循環不良；或皮膚知覺遲鈍者，需有人看護才可使用。
- 熱敷墊務必放入軟性外套內使用。
- 切勿在充滿高濃度氧氣之環境，以及接近有儲存或噴出氧氣之設備旁邊使用。
- 請勿強壓、硬折、扭擠本體。
- 請勿以電源線拉插頭或本體；或將電源線懸掛。
- 熱敷墊的本體部份有破裂或呈現發泡狀時，請勿再使用，並請更換新品。

TO REDUCE THE RISK OF BURNS, ELECTRIC SHOCK, AND FIRE, THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS.


- DO NOT USE WHILE SLEEPING.
- DO NOT USE ON AN INFANT.
- THIS PAD IS NOT TO BE USED BY OR ON AN INVALID, SLEEPING OR UNCONSCIOUS PERSON, OR A PERSON WITH POOR BLOOD CIRCULATION UNLESS CAREFULLY ATTENDED.
- DO NOT USE ON AREAS OF INSENSITIVE SKIN.
- BURNS CAN OCCUR REGARDLESS OF CONTROL SETTING, CHECK SKIN UNDER PAD FREQUENTLY.
- NEVER USE PAD WITHOUT COVER IN PLACE.
- DO NOT USE IN OXYGEN ATMOSPHERE.
- DO NOT USE PINS OR OTHER METALLIC MEANS TO FASTEN THIS PAD IN PLACE.
- DO NOT SIT ON OR CRUSH PAD AVOID SHARP FOLDS.
- NEVER PULL THIS PAD BY THE SUPPLY CORD.
- DO NOT USE THE CORD AS A HANDLE.
- CAREFULLY EXAMINE INNER COVER BEFORE EACH USE, DISCARD THE PAD IF INNER COVERING SHOWS ANY SIGN OF DETERIORATION.
- READ AND FOLLOW ALL INSTRUCTIONS ON BOX OR PACKED WITH PAD BEFORE USING.

使用方法


OPERATION

A. 濕式熱敷法：

FOR MOIST HEAT APPLICATION

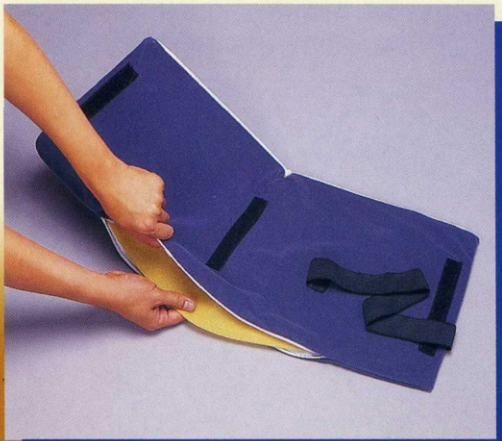

- 將海綿片浸水，多餘水份絞乾後，再將海綿片套在外套與熱敷墊本體中間。
- SATURATE SPONGE WITH WARM WATER, WRING OUT EXCESS, INSERT PAD & SPONGE INTO CLOTH COVER.

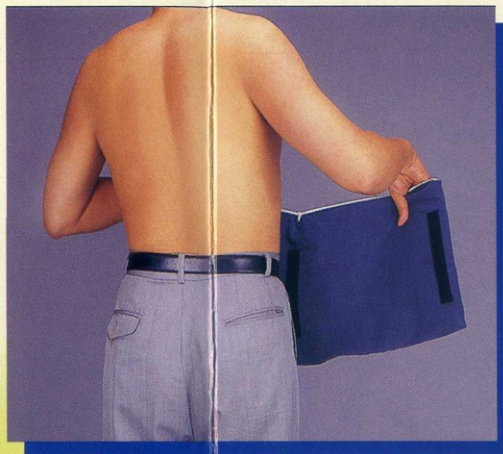


- 將海綿片的這一面，貼靠在身體部位上。
- PLACE THE SIDE OF THE PAD WITH THE SPONGE AGAINST BODY.

B. 乾式用法：

FOR DRY HEAT APPLICATION

- 用法同上（不需用到吸水海綿片）
- FOLLOW IN THE EXAMPLES ABOVE (YOU DON'T NEED TO USE THE ABSORBENT SPONGE)

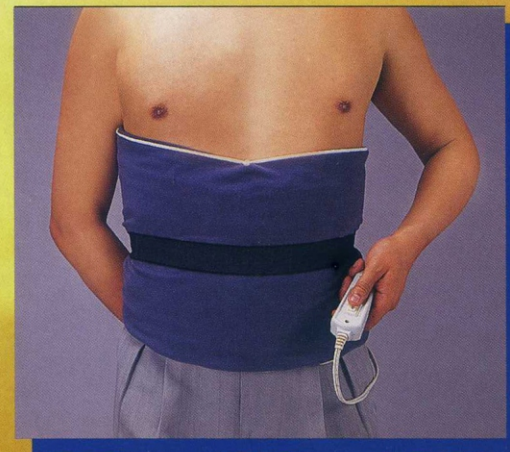


- 接入電源後，調整開關至適當溫度。
- PLUG INTO A POWER OUTLET & SET CONTROL TO DESIRED HEAT SETTING.

各種熱敷方式


HEATING WITH A VARIDE WRAPPING

■ 運用拉鍊、鬆緊帶、和黏扣帶的設計，可提供身體上各種不同部位之熱敷方式。


THE ZIPPER, ELASTIC BAND, HOOK & LOOP TAPES IS DESIGNED TO WRAP AROUND ANY PART OF BODY, PLACE HEAT PAD ON THE DESIRED AREA.

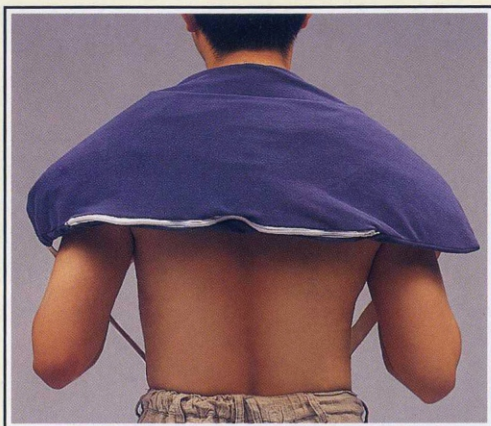


▲ 頸、肩膀 NECK, SHOULDERS

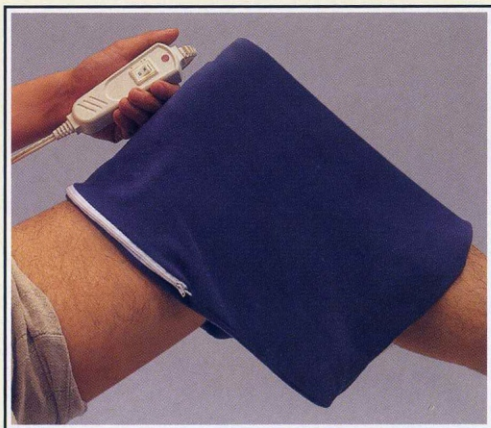

▲ 關節 JOINT

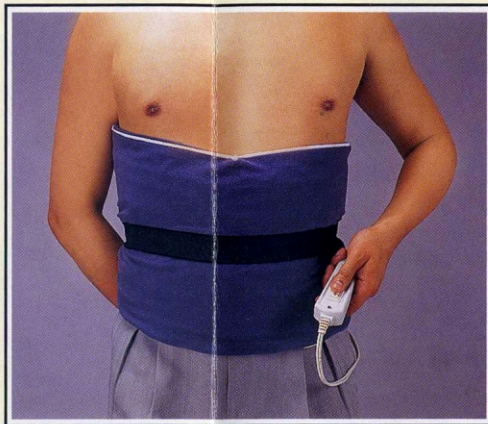

▲ 腰圍 LOWER BACK

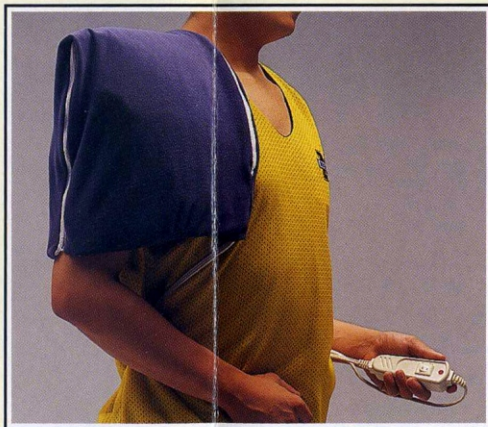


▲ 側肩 SIDE SHOULDERS

清洗外套說明：

■ 此產品外套可用冷水配合天然性洗衣劑來加以清洗，再風乾。

COVER WASHING INSTRUCTIONS

WASH CLOTH COVER IN COOL WATER ONLY, USING A MILD DETERGENT, USE COOL CYCLE TO DRY.

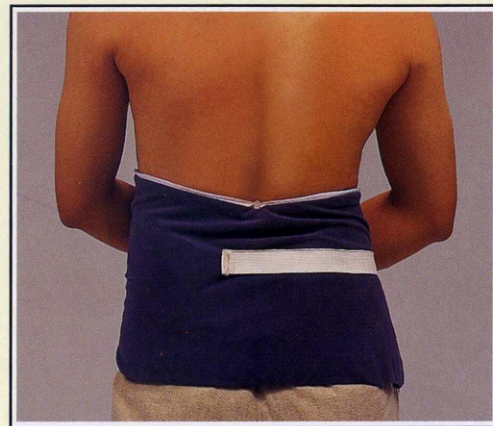

▲ 腹部 ABDOMEN



▲ 腳 FOOT

保證卡

憑本卡享受保證優待

購買日期	年 月 日	製造號碼
經銷店章		◆ 附記 1. 在正常使用下，如發生故障，自購買日起保證一年，修、換材料免費。 2. 因使用錯誤，導致故障，需換零件時，應繳付材料費。 3. 如故障請送原經銷店。 4. 本卡無購買日期或未蓋經銷店章者無效。

GUARANTEE

THIS PRODUCT IS WARRANTED FOR ONE YEAR FROM THE DATE OF PURCHASE AGAINST ELECTRICAL & MECHANICAL DEFECTS IN MATERIAL & WORKMANSHIP.

IF FOR ANY REASON THIS APPLIANCE SHOULD NOT OPERATE SATISFACTORILY, TEST THE SOCKET WHICH THIS APPLIANCE HAS BEEN CONNECTED TO WITH ANOTHER APPLIANCE. ALSO MAKE SURE THE PRONGS OF PLUG ON APPLIANCE ARE IN GOOD CONTACT WITH THE WALL OUTLET. IF YOU ARE SURE THE APPLIANCE IS DEFECTIVE, RETURN IT, TO APPLIANCE SERVICE STATION, FOR REPAIR OR REPLACEMENT WITHOUT CHARGE AT THE MANUFACTURER'S OPTION.

IF REPAIRS ARE NECESSARY AFTER GUARANTEE PERIOD A NOMINAL CHARGE WILL BE MADE.

DAMAGE TO APPLIANCE BY REASON OF ABUSE OR TAMPERING IS NOT COVERED BY THIS GUARANTEE. THIS GUARANTEE DOES NOT INCLUDE THE INSIDE OR OUTSIDE COVERING. WITH PROPER CARE THE COVERINGS SHOULD LAST INDEFINITELY.

中藍行實業有限公司
C-BLUE COMPANY LTD

高雄市三民區義華路 299 巷 8 號

NO. 8, LANE 299, YI HUA ROAD. 807 KAOHSIUNG, TAIWAN.

TEL: 886-7-3842438 FAX: 886-7-3840786